Bedienungsanleitung
OPERATING MANUAL
NOTICE D'EMPLOI
BEDIENINGSHANDLEIDING
MANUAL DE OPERACIÓN

Elektronische Säuglings- und Kleinkindwaage
Electronic baby & toddler scale
Pèse-bébé/bambin électronique
Báscula electrónica de bebés y niños pequeños
Elektronische weegschaal voor baby's & peuters

M112800

M112800-160316-Rev002-UM-de-e-es-fr-nl

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 www.ade-germany.de

Introduction

This symbol means: "Consult operating instructions". To achieve accurate results, read the instructions for use carefully before use and follow all instructions provided therein. Keep the instructions for use for future reference.

Your ADE electronic baby & toddler weighing has been manufactured for monitoring the weight of babies. The weighing result is used as an essential indication of the weight development of neonates and babies. Make sure the baby is only weighed following strictly the instructions of this Operating Manual and is never left unattended.

In order to achieve the best possible results it is imperative to read and follow the Operating Manual in all respects. Always keep the Manual in a safe place next to the scale.

Scope of delivery

- 1 x weighing scale base
- 1 x weighing tray
- 1 x weighing plate for toddlers
- 2 x AA batteries
- Operation manual

Basic safety precautions

- Handle your scale with care - it is a precision instrument.
- Never leave the infant / child unattended. Adult supervision is required at all times.
- Use the scale on a firm and level surface. In any case, make sure that the scale is secured against accidental dropping.
- Do not let the infant / child jump around on the scale.
- This product is not a toy and should not be used as one. When not in use, store out of reach of children.
- Before use, make sure that the mains voltage and type of current indicated on the rating plate match the mains voltage and type of current at the place of use.
- Make sure that the mains cable is not a potential trip hazard.
- Operate the unit only within the permissible ambient conditions.
- Do not operate the unit in hazardous or unstable environments.
- Locate the unit as far away as possible from equipment that generates electromagnetic or other interference.
- Portable and mobile RF communications equipment can affect electrical equipment. Before using high-frequency devices, observe their operating instructions.
- Use only approved accessories.
- To avoid the risk of fire and electric shock, use only ADE plug-in power supplies.
- Disconnect the mains plug before cleaning the appliance.
- Clean with a dry cloth, do not use harsh cleaning agents.
- Do not immerse the appliance in water or other liquids.

Safe operation is no longer guaranteed

- If the adapter shows visible signs of damage
- If the adapter no longer functions
- After extended storage in humid areas

In these cases, contact your service partner for your own safety!

Set up of the scale

The scale consists out of 3 parts – weighing base, baby weighing tray and toddler weighing plate.

Weighing scale base


Weighing tray

Weighing plate for toddler

The baby weighing tray or the toddler weighing plate is fixed and released via the turning knob underneath the baby weighing base. Please make sure, that the tray or plate is completely and securely fixed before use.

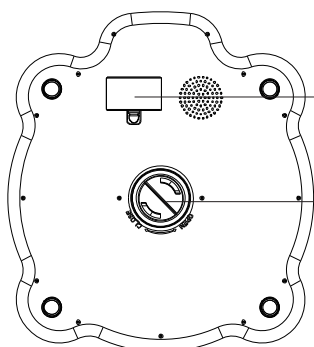Battery compartment

Lock knob

Power supply

The scale is powered by 2 x AA alkaline batteries. Turn the scale around and open the battery compartment. Insert the enclosed AA alkaline batteries. Please observe the correct polarity as shown inside the battery compartment.

Replace the batteries when „Lo“ is being displayed. It indicates low battery power or use the mains connection.

Control and Display

Control

Keys	Description
ON / OFF	Turn on and zero setting – automatic switch off
HOLD	Activates the auto hold function – to unlock press again
TARE	Activates the tare function

Display

Display symbol	Description
◀	Zero position
Hold	Auto hold
kg lb	Weighting units kg or lbs

Operation

Weighing

1. Start the scale by pressing the **ON/OFF** button. After a function test the last measured weight will be displayed for 3 seconds followed by 0.00 and ◀. The scale is ready for use and is set automatically to zero.
2. The baby can now be placed on the scale and the weight will be shown on the display in 10 g steps.

Auto hold function

Place the baby on the scale. The display will flush until the actual weight has been taken and **HOLD** will be shown on the display.

To cancel the hold function press the **HOLD** button on the key board.

The scale will switch of after approx. 60 seconds, when not in use.

Tare function

The tare function can be used to disregard an additional weight e.g. a nappy or towel on the weighing tray surface. The tare function is operated as follow: Place the object on the tray surface and press the **TARE** button. The additional weight will be subtracted and 0.00 will be shown on the display. A baby can now be placed on the scale surface and the net weight will appear on the display. The maximum tare weight must not exceed 500g.

Unit switch function

Press the **ON/OFF** button twice during the function test. By pressing again the **ON/OFF** button you can select the requested units **kg** or **lbs**. The **kg** or **lb OZ** symbol will appear in the display to indicate the chosen weighing unit.

Recall of memorized weight

Start the scale via the **ON/OFF** button and the last measured weight will be displayed for approx. 3 sec. in order to determined weight changes of the infant.

Cleaning & maintenance

Cleaning

Clean your scale with a soft moist towel or conventional disinfectant only. Do not use any aggressive cleaning solutions.

Fault messages

1. Why does the display read 'lo'?

Low battery. Replace the batteries.

2. The message displays 0-Ld

Overload warning; remove the weight immediately; otherwise, permanent damage to the scale will occur.

Storage and transport

Keep all original parts and pieces for an eventual sending back of the scale in order to avoid any possible damage during transport which automatically leads to no warranty.

Disconnect all cables before transport to avoid unnecessary damages.

Waste disposal

Electronic devices are not domestic waste.

If the scale is out of use the responsible waste disposal association will be glad to advise you on the measures necessary to dispose of it properly.

Do not throw used batteries into domestic waste. Please inform yourselves at the responsible waste disposal association how to dispose used batteries.

Warranty

A two year warranty period from date of purchase will be guaranteed for either repair or replacement once it is proven that the defects are caused by poor material or poor workmanship (please keep the receipt for presentation). All moveable parts such as batteries, cables, plug-in power supply units, rechargeable batteries etc. are excluded. Warranty does not apply to normal appearance of usage or sign of wear or any damages caused by accidents or misuse. Eventual warranty- or liability claims are subject to using the original ADE accessories or spare parts. Products which have been opened by unauthorized persons are automatically out of warranty. Customers abroad are requested to contact their local dealers in the event of any warranty becoming applicable.

Technical data

Capacity:	50 kg
Graduation:	10 g
Weighing unit selectable:	kg, lb
Battery operation	2pcs 1.5V AA alkaline batteries

Weight & Dimensions:

Scale (LxWxH):	631 x 362 x 111 mm
Baby tray (LxWxH):	631 x 320 x 60 mm
Weight:	approx. 3,1 kg

Declaration of conformity by the manufacturer

These scales have been manufactured in accordance with the harmonised European standards. They conform to the regulations of the EC-directives listed below:

- EMC Directive 2014/30/EU
 - Low Voltage Directive 2014/35/EU
 - RoHS Directive 2011/65/EC
- in the currently valid versions.

This declaration loses its validity if modifications are made to the scale without our approval.

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ADE Germany GmbH